

Technique Creates Beauty

Dual-Max

IPL SHR Hair Removal Skin Treatment Machine

F
A
S
T

S
E
R
V
I
C
E

A
N
D

B
E
S
T

V
A
L
U
E





Dual-Max

FDA-Cleared IPL Machine
IPL SHR Hair Removal Skin Treatment Machine

■ Application

- Permanent Hair Removal
- Photo Rejuvenation
- Pigmentation Correction
- Vascular Therapy
- Acne Treatment
- Wrinkle Removal

Standard



E-light:

Spot size: 15*35mm²
Wavelength: 430/530/640-1200nm
Function: Photo Rejuvenation, Pigmentation Correction, Vascular Therapy, Acne Treatment, Wrinkle Removal



SHR:

Spot size: 15*50mm²
Wavelength: 650-950nm for SHR
Function: Fast Hair Removal, suitable for skin types I to V

Optional



IPL:

Spot size: 15*50mm²
Wavelength: 430/530/640-1200nm
Function: Photo Rejuvenation, Pigmentation Correction, Vascular Therapy, Acne Treatment, Wrinkle Removal



IPL:

Spot size: 12*33mm²
Function: Photo Rejuvenation, Pigmentation Correction, Vascular Therapy, Acne Treatment, Wrinkle Removal



RF:

Spot size: 3 Tips
Function: Skin Lift & Tighten for Eyes, Face, Body




Technique Creates Beauty

■ IPL SHR Work Theory



The IPL SHR unique gradual heating method effectively damages the hair follicles while avoiding injury to the surrounding skin. The sweeping In-motion delivery technique keeps the IPL in constant motion ensuring full coverage.

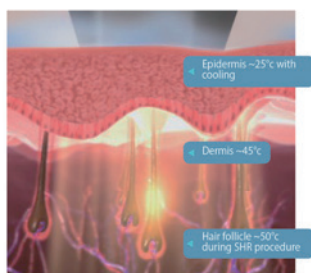

Simultaneous contact cooling prevents burns and allows for more comfortable treatment and a specialized applicator allows you to treat large areas in less time.

Whether you have light, dark or tanned skin, coarse or fine hair, SHR hair removal solutions achieve maximum results for all areas of the body with treatments that are safe, quick and best of all-virtually painless.

■ Proven Clinical Results



Hair Removal



Hair Removal



Acne Removal



Pigmentation Removal



Vascular Removal



Skin Rejuvenation

Technique Creates Beauty

■ Advantages



Quality Assurance

- ▶ TUV Medical CE, CFDA and American FDA-cleared, ensure safe operating and excellent quality.

Integrated Self-cooled FDP Handpiece

- ▶ Self-cooled system - pain free; Fast Delivery Pulse technology, up to 10Hz - shorten operating time; Different spot size - work on face or body conveniently.

Germany Imported Xenon Lamp

- ▶ Germany Heraeus Xenon Lamp lifetime 800,000 flashes, stable energy output and longer life span.

Japan Imported Capacitors

- ▶ High performance Japan Nippon Capacitors 12000uF x 3 , ensure equal energy output of each pulse at the speed of 10 shots per second.

Intelligent Software

- ▶ Simple mode for beginners, can get Pre-set parameters automatically according to skin types, hair density, colors and thickness.
- ▶ Professional mode for experts, can set different parameters according to your patients.

Technique Creates Beauty

Why Choose Dual-Max for Your Next IPL (Photofacial) Treatment?

Dual-Max IPL is the gold standard for photorejuvenation. This FDA-cleared, TUV Medical CE-cleared, unique technology provides dramatic results for a variety of benign conditions, including age spots, sun-induced freckles, symptoms of rosacea, unsightly veins, and other blemishes.

What Conditions does IPL Photorejuvenation treat?

- Young patients with milder degrees of sun damage or early aging changes
- Rejuvenation of skin of the face, neck, chest, shoulders, arms or hands
- Pigment Changes/Sun Damage or Aging Skin
- Ruddy Skin Appearances ("Dyschromia") - rosacea, redness of the chest and neck
- Enlarged Pores
- Fine Lines
- Telangiectasia (broken capillaries)
- Unwanted Freckles
- Hyperpigmentation or brown spots

What to do Before your IPL Treatment:

- Avoid tanning and excessive sun exposure 4 weeks prior to your treatment.
- The treatment areas should be cleansed and all makeup or creams removed.
- Each treatment takes approximately 15 to 45 minutes depending on the size of the area treated.

What to do After your IPL Treatment:

- A sunburn sensation is common. Cold packs may be used for ten minutes after your treatment.
- The skin may appear pink and flushed. This usually resolves after several hours.
- Age spots and melasma patches will darken over several days and either exfoliate off over 7 to 10 days or fade over several weeks.
- Striping is common in areas treated off the face, like the arms and chest, and require follow up treatment to even out the results.
- Gentle cleansing and sun protection are important while the skin is sensitive.
- It is very important to avoid tanning 4 weeks before and 4 to 6 weeks after your treatments.
- In general, there is little downtime with IPL. Occasional bruising and swelling may occur but is temporary and resolves over 3 to 7 days. You can resume normal activities immediately after your treatment.
- Gradual improvement will be noted over your course of treatments.

Technique Creates Beauty

■ IPL Treatment Menu Pricing

IPL Photofacials Pricing

Full Face: \$200

Face and Chest: \$300

Face, Neck, and Chest: \$380

Hands: \$140

IPL Pigmentation Pricing

Such as redness, freckles and post-trauma pigmentation such as acne scars and age spots.

Nose: \$40

Cheeks: \$70

Quarter Face: \$83

Half Face: \$111

Full Face: \$139

Neck: \$98

Decolletage (Chest): \$139

Face + Neck: \$213

Face + Neck + Decolletage: \$279

Hands: \$139

Shoulders: \$139

Lower Arms: \$174

Full Arms: \$347

Lower Legs: \$272

Full Legs: \$544

Upper Back: \$376

Full Back: \$531

IPL Hair Removal Pricing (per treatment)

Lip: \$59

Chin: \$69

Jaw: \$79

Sideburns: \$99

Full Face: \$250

Ears: \$49

Brow Middle: \$30

Brow Forehead: \$49

Between Breast: \$49

Areola: \$59

Lower Navel Below Belt Line + Belly Line: \$79

Bikini Line Traditional 2 Fingers From Crease: \$139

Bikini line Extended with Inner Thighs: \$169

Brazilian Everything Front + Back: \$225

Buttocks (Glutes): \$69

Brazilian Backside: \$69

Inner Thighs: \$69

Underarms: \$119

Upper Arms: \$129

Lower Arms: \$129

Full Arms: \$199

Hands: \$99

Feet: \$69

Back of Neck: \$69

Lower Back: \$69

Upper Legs: \$219

Lower Legs: \$219

Full Legs: \$399

Packages of 8 Receive 50% Off

Technique Creates Beauty

■ Technique Specification

Model NO.

IPL Peak Power	3500W
Wavelength	SHR 650nm-950nm; Elight 430/530/640-1200nm standard; 480/560/590/690/750-1200nm optional
Energy Density (Fluence)	1-26J/cm ² (SHR), 1-60J/cm ² (Elight)
Spot Size	15×50mm ² (SHR), 15×35mm ² (Elight) (12×33mm, 10×40mm optional)
Pulse Duration	1-15ms
Pulse	Single
Pulse Repetition Rate	1-10Hz (SHR)
Timer	1-30ms
Cooling	Air + Water + Sapphire contact cooling (-4℃~+4℃)
Stand-by Working	Continuously for 12 hours
Display	8.4" True Color LCD Touch Screen
Net Weight	46KG
Dimensions (W x L x H)	54×50×118cm

